# Lab 3: Reading and Manipulating Short Reads

# Martin Morgan\*

# 14-18 June, 2010

# Contents

1	1 Introduction	1
2	2 Aligned read input	2
	2.1 Navigating Solexa output	 . 2
	2.2 readAligned and the AlignedRead class	 . 3
	2.3 Subsets and filters	 . 7
	2.4 Cautions	 . 9
3	3 Additional input functions	9
4	4 Quality assessment	11
	4.1 Generating a QA report	 . 11
	4.2 Exploring qa	
	4.3 Frequent sequences	 . 13
	4.4 Cycle-specific qualities and base calls	 . 16
5	5 Coverage	18
6	6 Summary	18
7	7 Session information	19

# 1 Introduction

This practical uses the *ShortRead* package to input aligned and other short read data files, and illustrates some of the available sequence manipulation and quality assessment tools. Activities during the lab are posed as exercises. So, as a first exercise:

## Exercise 1

Start an R session, and load the ShortRead package.

- > library(ShortRead)
- > packageDescription("ShortRead")\$Version

 $<sup>^*\</sup>mathrm{Fred}$  Hutchinson Cancer Research Center, Seattle, WA 98008

```
[1] "1.6.2"
```

Confirm that the version of your package is at least as recent as the version in this document. Seek assistance from one of the course assistants if you need help getting the current version of ShortRead.

The course also requires access to sample data.

#### Exercise 2

Obtain and load the CSAMA10 package. Use

> install.packages("path/to/CSAMA10.tar.gz", repos=NULL)

# 2 Aligned read input

This section illustrates input of aligned reads. It focuses on aligned reads produced by the Solexa Genome Analyzer ELAND software; reading data produced by software such as MAQ or Bowtie, or in BAM format is described in the ShortRead 'Overview' vignette and on the readAligned help page.

# 2.1 Navigating Solexa output

Vendor and third-party software is likely to process raw images, base calling (*Rolexa* is a *Bioconductor* package providing alternative base calling), and alignment to a reference genome (see *BSgenome* and the matchPDict function of *Biostrings*). *Bioconductor* packages might enter a typical work flow after alignment. The following starts with reads aligned with *ELAND*, an alignment program from the Illumina Genome Analyzer II (GAII) platform.

We'll start by creating a variable, extdataDir, containing the directory holding the sample GAII data, and check to make sure that we have defined the right path.

A key functionality provided by *ShortRead* is input of a diversity of file types, both of files from the Illumina pipeline and from other software and in formats appropriate for other technologies. The interface to these input functions is meant to facilitate reading one or more files into a single object. The interface is like that for list.files: provide a directory path where relevant files are to be found, plus a regular expression to select the files you are interested in.

### Exercise 3

Use list.files to display all files in the extdataDir directory.

```
> list.files(extdataDir)
```

We'll select just one file to use in this analysis, based on files matching the regular expression .\*\_export\_head.txt. These files are produced using ELAND software run in eland-extended mode. This mode produces files, one for each lane (and 'end' of paired end reads) that summarize diverse features of all reads, and is a very convenient starting point for analysis; other inputs (e.g., MAQ, Bowtie, BWA, or BAM files) are also supported.

#### Exercise 4

We'll use an abbreviated file for most parts of this lab. The abbreviated file contains the first 500,000 reads from a single lane of an Illumina paired-end run. It is from lane 1, and we will use the first end only. The name of the file is  $s_1_1$ export\_head.txt. We will use this as the 'pattern' to match, and check that we have specified the pattern appropriately

```
> pattern <- "s_1_1_export_head.txt"
> list.files(extdataDir, pattern)
[1] "s_1_1_export_head.txt"
```

Our success shouldn't be too surprising in this case, but it often pays to check. For instance, the pattern above would also match a file with the same name but with .tar.gz appended!

## 2.2 readAligned and the AlignedRead class

The readAligned function can be used to input aligned reads. The first argument is a directory path where alignment files are to be found. The second argument is the regular expression to select files to be read. An optional third argument allows the user to specify which type of file is to be read in.

### Exercise 5

Use readAligned to read in our abbreviated version of lane 1, specifying Solex-aExport as the type of file.

```
> aln <- readAligned(extdataDir, pattern, "SolexaExport")</pre>
```

See the help page for readAligned for additional details about supported file types.

What does readAligned input? It inputs the short read sequences and base call qualities, and the chromosome, position, and strand information associated with short read alignments. This information is expected to be provided by all short read alignment software.

#### Exercise 6

Display the object we have read in.

> aln

```
class: AlignedRead
```

length: 500000 reads; width: 35 cycles

chromosome: 255:255:255 255:255:255 ... QC QC

position: NA NA ... NA NA strand: NA NA ... NA NA alignQuality: NumericQuality

alignData varLabels: run lane ... filtering contig

There are 500000 reads in the object, each read consisting of 35 nucleotides. View the first several reads and query information about, e.g., the number of reads that align to each strand, or the number of positions recorded as NA.

### > head(sread(aln))

```
A DNAStringSet instance of length 6 width seq
```

- [1] 35 AAAAAAAAAAAAAAAAAAAAAAAAAAAAAAAAA
- [2] 35 AAAAAAAAAAAAAAAAAAAAAAAAAAAAAAAA
- [3] 35 AAAAAAAAAAAAAAAAAAAAAAAAAAAAAAAAA
- [4] 35 AAAAAAAAAAAAAAAAAAAAAAAAAAAAAAA
- [5] 35 AAAAAAAAAAAAAAAAAAAAAAAAAAAAAAAA
- [6] 35 GCAAGTTAAGAAGAGAGCAGAAGAACGTTTTTA

```
> table(strand(aln), useNA="ifany")
```

> sum(is.na(position(aln)))

### [1] 261002

Notice how elements of the alm object are extracted using accessors such as sread and strand; these are described on the help page for the class of the alm object (indicated in the display of alm, above, as class AlignedRead); note that the help page refers to the help page for accessors to enumerate additional ways of accessing the data.

What are all the NA values returned by strand and position? These correspond to reads that did not align to the reference genome used by *ELAND*; that about 1/2 the reads do not align is below normal, making this an interesting opportunity for quality assessment. The strand function returns a factor with three levels. The first two describe reads aligned to the plus and minus strands, the third (\*) is available for successful alignments where strand information is irrelevant.

Aligned reads contain several different kinds of information about 'quality'. Individual bases are assessed for quality during base calling. These 'raw' base qualities are 'calibrated' during *ELAND* alignment; details of calibration are to be found in Illumina documentation. The alignments themselves also have qualities associated with them, with the details of alignment quality differing between alignment algorithms.

#### Exercise 7

Retrieve calibrated base quality from aln.

```
> head(quality(aln))
class: SFastqQuality
quality:
  A BStringSet instance of length 6
    width seq
Γ17
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
Γ27
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
[31]
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
[4]
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
[5]
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
Γ67
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
```

These qualities are string-encoded  $-10\log_{10}$  probabilities. The encoding in this case follows a convention established by Illumina. The details of the encoding can be obtained by querying quality(aln) for its alphabet; the letter A corresponds to a  $-10\log_{10}$  score of 1.

Numeric values are readily retrieved as a matrix, with rows corresponding to reads and columns to cycles. Computations can then be performed on them, e.g., to determine average calibrated quality scores as a function of cycle.

```
> alf <- alphabet(quality(aln))
> m <- as(quality(aln), "matrix")
> colMeans(m)

[1] 21.45995 19.91676 17.98661 19.67178 18.27315 18.61814
[7] 18.99522 18.81180 18.69941 18.47509 19.12563 17.39689
[13] 18.26270 18.85640 17.82326 17.43106 17.95288 17.01137
[19] 17.56018 17.73439 16.51199 16.99086 17.80355 16.38983
[25] 17.35305 17.15449 15.99980 16.38890 16.59765 16.18294
[31] 13.48905 13.03266 13.00785 12.75729 11.96461
```

Alignment qualities are accessible with alignQuality. This returns an object that can contain quality scores in different formats; to extract the actual quality scores, use quality. Reads failing to align or to align in multiple locations have an alignment quality of 0.

#### Exercise 8

Retrieve the alignment quality scores, determine how many align poorly, and visualize the distribution (Figure 1) of scores.

```
> alignQuality(aln)
class: NumericQuality
quality: 0 0 ... 0 0 (500000 total)
> q <- quality(alignQuality(aln))
> sum(q==0)
[1] 295387
```

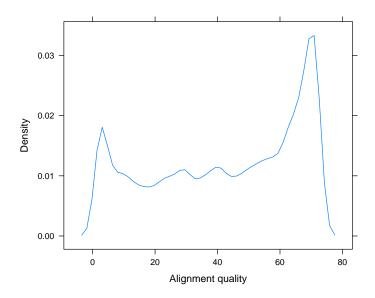


Figure 1: Alignment quality

```
> print(densityplot(q[q>1], plot.points=FALSE,
+ xlab="Alignment quality"))
```

Alignment algorithms produce information in addition to basic data about chromosome, position, and strand alignment. The exact content varies between algorithms, and is available with alignData. alignData returns an Aligned-DataFrame object that contains these data and a metadata description of them. For instance, ELAND includes information about whether the read passed a base-calling filter (based on strength and consistency of early bases in the read), in addition to the lane, tile, x and y coordinate of each read.

### Exercise 9

Use the alignData function to extract the additional information in the ELAND alignment file. The underlying data in this object can be accessed as though it were a data frame, for instance to tally the number of reads passing Illumina base calling filter.

```
> alignData(aln)
```

```
An object of class "AlignedDataFrame"
readName: 1, 2, ..., 500000 (500000 total)
varLabels and varMetadata description:
run: Analysis pipeline run
lane: Flow cell lane
...: ...
contig: Contig
(7 total)
```

```
> table(alignData(aln)$filtering)
```

Y N 287222 212778

### 2.3 Subsets and filters

A very common operation is to reduce the number of reads used for subsequent analysis. This can be done in a coordinated fashion by creating a subset of aln.

### Exercise 10

Select just the aligned reads passing Illumina filtering, and aligning to the reference genome.

```
> filtIdx <- alignData(aln)$filtering=="Y"
> alignedIdx <- !is.na(strand(aln))
> aln[filtIdx & alignedIdx]

class: AlignedRead
length: 197432 reads; width: 35 cycles
chromosome: chr11.fa chr9.fa ... chr8.fa chr4.fa
position: 104853312 3036336 ... 44295163 47191474
strand: - - ... - -
alignQuality: NumericQuality
alignData varLabels: run lane ... filtering contig
```

A different approach to subsetting is to use objects of class SRFilter. These can be particularly useful as an argument to readAligned, in addition to use in interactive sessions.

### Exercise 11

Construct instances of built-in filters to select reads passing the Illumina filtering criterion, and uniquely aligning to a fully assembled chromosomes. These can be 'composed' into a single overall filter, and applied to restrict available reads.

```
> filt1 <- alignDataFilter(expression(filtering=="Y"))
> filt2 <- chromosomeFilter("chr[0-9XYM]+.fa")
> filt <- compose(filt1, filt2)
> caln <- aln[filt(aln)]
> caln

class: AlignedRead
length: 195719 reads; width: 35 cycles
chromosome: chr11.fa chr9.fa ... chr8.fa chr4.fa
position: 104853312 3036336 ... 44295163 47191474
strand: - - ... -
alignQuality: NumericQuality
alignData varLabels: run lane ... filtering contig
```

The filters developed above could be used to filter reads while being read in to R, e.g,. with

```
> readAligned(extdataDir, pattern, type="SolexaExport", filter=filt)
```

The srFilter function can be used to create custom filters. The idea is that filter functions accept a single argument x that is an object to be filtered, and returns a logical vector that can be used to select elements of the object.

#### Exercise 12

As a first example, write and use a filter to select only a single read from all that align to a particular chromosome and position (normally one would also like to include strand as a uniqueness criterion; see the occurrenceFilter function for a more complete implementation).

```
> ualignFilter <- srFilter(function(x) {
      ## create a numerical index of reads. Divide the index, and
      ## position information between chromosomes. Select the index of a
      ## single read at each unique position. Return the selected index
      ## as a logical vector with the same length as x
      oindex <- seq_len(length(x))</pre>
      index <- tapply(oindex, chromosome(x), c)</pre>
      pdup <- tapply(position(x), chromosome(x), duplicated)</pre>
      keep <- oindex %in% unlist(mapply(function(i, p) {</pre>
          i[!p]
      }, index, pdup))
+ }, name="select only one read per position")
> caln[ualignFilter(caln)]
class: AlignedRead
length: 188219 reads; width: 35 cycles
chromosome: chr11.fa chr9.fa ... chr8.fa chr4.fa
position: 104853312 3036336 ... 44295163 47191474
strand: - - ... - -
alignQuality: NumericQuality
alignData varLabels: run lane ... filtering contig
```

The filter functions built-in to *ShortRead* use a 'factory' pattern to create instances of each filter that 'remember' how the filters were created. For instance, chromosomeFilter("chr2.fa") creates an instance of the chromosome filter to select only chromosomes matching chr2.fa.

#### Exercise 13

As an advanced example, the following filter subsamples a (user-specified) number of reads. The samplingFilter function uses the factory pattern, so filters created with it remember how many reads to sample.

```
> samplingFilter <- function(sampleSize) {
+     srFilter(function(x) {
+         idx <- seq_len(length(x))
+         idx %in% sample(idx, sampleSize)
+     }, name="Demo sampling filter")
+ }
> sample100 <- samplingFilter(100)
> caln[sample100(caln)]
```

class: AlignedRead

length: 100 reads; width: 35 cycles

chromosome: chr16.fa chr4.fa ... chr17.fa chr13.fa position: 39017877 83241342 ... 47359791 106761505

strand: - + ... + +

alignQuality: NumericQuality

alignData varLabels: run lane ... filtering contig

### 2.4 Cautions

There are several confusing areas associated with reading data aligned with various software packages. (1) Some alignment programs and genome resources start numbering nucleotides of the subject sequence at 0, whereas others start at 1. (2) Some alignment programs report matches on the minus strand in terms of the 'left-most' position of the read (i.e., the location of the 3' end of the aligned read), whereas other report 'five-prime' matches (i.e., in terms of the 5' end of the read), regardless of whether the alignment is on the plus or minus strand. (3) Some alignment programs reverse complement the sequence of reads aligned to the minus strand. (4) Base qualities are sometimes encoded as character strings, but the encoding differs between 'fastq' and 'solexa fastq'. It seems that all combinations of these choices are common 'in the wild'.

The help page for readAligned attempts to be explicit about how reads are formatted. Briefly:

- Subject sequence nucleotides are numbered starting at 1, rather than zero. readAligned adjusts the coordinate system of input reads if necessary (e.g., when reading MAQ alignments).
- ELAND and Bowtie alignments on the minus strand are reported in 'left-most' coordinates systems.
- ELAND and Bowtie alignments on the minus strand are not reverse complemented.
- Character-encoded base quality scores are interpreted as the default for the software package whose output is being parsed, e.g., as 'Solexa fastq' for *ELAND*. The object returned by quality applied to an *AlignedRead* object is either *FastqQuality* or *SFastqQuality*.

Alignment programs sometimes offer the opportunity to customize output; such customization needs to be accommodated when reads are input using *Short-Read*.

# 3 Additional input functions

Many sequence alignment programs are producing 'bam' (Binary AlignMent) files. These can be read in to R using readAligned with type="BAM". readAligned filters reads to include only those reads alighed without indels; more flexible representation is available with readGappedAlignments, which inputs the reads and their 'cigar' representation. A cigar is a short run length encoding providing information about how the sequence of nucleotides in the read align to the

reference, including information about indels. Even greater flexiblity is provided by scanBam function in the *Rsamtools* package. This will parse *BAM* files that are either local or remote (e.g., from the 1000 genomes project), and can be used to efficiently select reads satisfying particular criteria, e.g., overlapping a particular range or aligned to the minus strand. See the help page for scanBam and ScanBamParam for additional information.

ShortRead, Biostrings, and the standard input functions from R provide additional tools for reading Illumina and other alignment formats. The read-XStringColumns function provides a convenient way to read DNA and quality sequences into compact data structures. readFasta and its counterpart readFasta provide tools for reading FASTA- or FASTQ- (i.e., including quality annotation) formatted files.

#### Exercise 14

Files \_sequence.txt contain fastq-formatted sequence and quality scores. A sample of this file type is available. Read these in to data structured defined in ShortRead. The content of reads can be retrieved with the accessor functions id, sread, and quality.

```
> reads <- readFastq(extdataDir, "s_1_1_sequence_head.txt$")</pre>
> head(id(reads))
  A BStringSet instance of length 6
    width seq
[1]
       24 HWI-EAS88_3:1:1:33:484/1
[2]
       24 HWI-EAS88_3:1:1:33:272/1
[3]
       24 HWI-EAS88_3:1:1:31:594/1
[4]
       24 HWI-EAS88_3:1:1:33:383/1
Γ51
       24 HWI-EAS88_3:1:1:35:216/1
[61
       25 HWI-EAS88_3:1:1:883:458/1
```

The readPrb functions reads 'raw' base call quality scores from \_prb files in the baseCallPath directory; the result is a *BStringSet* object that compactly represents the quality scores in a way analogous to the results of quality applied to aln.

#### Exercise 15

The \_export.txt files read by readAligned are tab-delimited text files. The goal of this exercise is to read the DNA sequence and quality score columns in to R as DNAStringSet and BStringSet objects. The DNAStringSet and BStringSet classes represent DNA or 'biological' strings; they extend the base class XString.

Start by parsing the first line of a \_prb file to get a sense of its content. Specify the colClasses to be imported, using NULL to indicate that a column should be skipped, and DNAString or BString to indicate columns that are to be read as the corresponding data types (classes). Read the file with readXStringColumns.

```
> fl <- list.files(extdataDir, pattern, full=TRUE)
> cols <- strsplit(readLines(fl, 1), "\t")[[1]]
> length(cols)
[1] 22
```

```
> cols[9:10]
[1] "AAAAAAAAAAAAAAAAAAAAAAAAAAAAAAA
[2] "ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU"
> colClasses <- rep(list(NULL), 22)
> colClasses[9:10] <- c("DNAString", "BString")</pre>
> strings <-
      readXStringColumns(extdataDir, pattern, colClasses=colClasses)
> head(strings[[2]])
  A BStringSet instance of length 6
    width seq
[1]
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
Γ21
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
[31]
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
[4]
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
Γ51
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
[61
       35 ZZZZZZZZZZZZZZZZZZZZZZZZZZZZZZUUUUU
```

For each column in colClasses, readXStringColumns parses the corresponding column in all files matching the function argument pattern into a single XStringSet, i.e., concatenating the content of all files into a single object.

# 4 Quality assessment

This part of the course addresses *ShortRead* facilities for assessing quality, primarily of Solexa data. The *ShortRead* functionality is mean to complement rather than replace QA tools provided by the *ELAND* pipeline.

# 4.1 Generating a QA report

Creating a QA report is a two-step process. The first step is to visit necessary files to collate information in a compact representation. The second step is to present the information in a useful format.

The qa function collates information for the QA report. It visits each \_export.txt file, and extracts information on reads and their qualities. Evaluation of the function is straight-forward, e.g., qa <- qa(extdataDir, ".\*\_export.txt", type="SolexaExport"). The process of collating files can be time consuming (each export file must be parsed, taking 3-4 minutes per file) and memory intensive (lanes are processed independently of one another, but processing a full lane consumes 2-3 GB of memory). The return value of the qa function is actually quite compact, and easy to work with.

# Exercise 16

Rather than collating information during the lab, we load the data from a previously stored instance.

```
> data("qa_080828_081110", package="CSAMA10")
> qa
```

```
class: SolexaExportQA(9)

QA elements (access with qa[["elt"]]):
  readCounts: data.frame(8 3)
  baseCalls: data.frame(8 5)
  readQualityScore: data.frame(12288 4)
  baseQuality: data.frame(752 3)
  alignQuality: data.frame(629 3)
  frequentSequences: data.frame(1200 4)
  sequenceDistribution: data.frame(4540 4)
  perCycle: list(2)
   baseCall: data.frame(1400 4)
   quality: data.frame(6664 5)
  perTile: list(2)
   readCounts: data.frame(7200 4)
  medianReadQualityScore: data.frame(7200 4)
```

One feature of *ShortRead* that can speed up this stage of the operation is the use of clustered computer resources and the *Rmpi* or *mullticore* package; qa uses the srapply function to automatically detect and distribute collation tasks across pre-established nodes. This is outlined in more detail in a subsequent section.

The QA information collated from the \_export.txt files is summarized into an html report using the report function. The command to create the report is rpt <- report(qa, dest=tempfile()). This creates an HTML file at the location specified by the argument dest.

### Exercise 17

Create a report from the qa object loaded in the previous exercise. Do this with the command

```
> rpt <- report(qa, dest=tempfile())
```

View the report in your browser with

```
> browseURL(rpt)
```

The QA report provides summary statistics about the numbers of reads and alignments, base calls and qualities, characteristics of per-lane and per-tile read quality, and other information. The QA report is self-documenting, providing a narrative description of each section.

# 4.2 Exploring qa

The qa object is a list-like structure with several entities.

# Exercise 18

The readCounts element of qa is a simple data frame summarizing, on a per-lane bases, the total number of reads, the number of reads passing Solexa internal filtering, and the number of aligned reads.

```
> qa[["readCounts"]]
```

```
read filtered aligned s_1_1_export.txt 3668433 2278945 1910666 s_2_1_export.txt 4230424 3239956 2771169 s_3_1_export.txt 4003465 3089375 1720396 s_4_1_export.txt 4521919 3446177 2571235 s_6_1_export.txt 4004807 3127297 1985855 s_7_1_export.txt 3546869 2732974 1368590 s_8_1_export.txt 4232977 3291627 2379709 s_5_1_export.txt 2842633 2399387 2320140
```

Lanes 1-4 and 6-8 correspond to biologically interesting samples; lane 5 is the Solexa  $\varphi X$ -174 control lane. The number of reads (between 2.8 and 4.5 million) is low for a typical experiment (official guidelines are provided in Solexa documentation).

It can be difficult to scan large numbers, so the QA report template defines functions that help to display the information in a more comprehensible fashion. Source these files into your current R session, and read the second and third columns as a proportion of the first.

> ShortRead:::.ppnCount(qa[["readCounts"]])

```
read filtered aligned
s_1_1_export.txt 3668433
                            0.621
                                     0.521
                                     0.655
s_2_1_export.txt 4230424
                            0.766
                                     0.430
s_3_1_export.txt 4003465
                            0.772
s_4_1_export.txt 4521919
                            0.762
                                     0.569
s_6_1_export.txt 4004807
                            0.781
                                     0.496
s_7_1_export.txt 3546869
                            0.771
                                     0.386
s_8_1_export.txt 4232977
                             0.778
                                     0.562
s_5_1_export.txt 2842633
                            0.844
                                     0.816
```

Refer to the QA report for further commentary on this and other aspects of the report

### 4.3 Frequent sequences

A feature of raw reads, and of many subsequent stages of short read analysis, is a power law-like relationship between the number of times a read occurs, and the number of occurrences of a particular sequence in the sample. This information is contained in the sequenceDistribution element of qa, and is the result of running the tables command (defined in ShortRead) on a DNAStringSet object.

# Exercise 19

Retrieve the sequenceDistribution element from qa; it is a simple data frame. Look at the contents of the data frame using head, and select just lane 5 raw reads. plot the power-law relationship between the number of reads and number of times reads occur. As an alternative display of the same information, plot the cumulative number of reads as a function of the number of times a read occurs.

```
> df <- qa[["sequenceDistribution"]]
> df5raw <- df[df$lane=="s_5_1_export.txt" & df$type=="read",]
> head(df5raw)
```

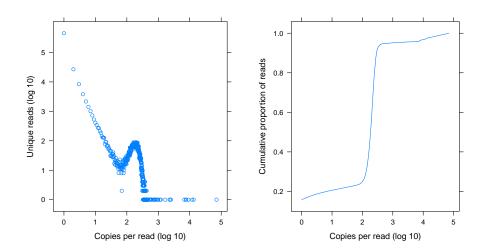


Figure 2: Number of copies of unique reads

```
nOccurrences nReads type
                                            lane
3337
                1 450124 read s_5_1_export.txt
3338
                   26810 read s_5_1_export.txt
3339
                3
                    8366 read s_5_1_export.txt
3340
                4
                    3806 read s_5_1_export.txt
3341
                5
                    2166 read s_5_1_export.txt
3342
                6
                     1417 read s_5_1_export.txt
 print(xyplot(log10(nReads)~log10(nOccurrences), df5raw,
               xlab="Copies per read (log 10)",
               ylab="Unique reads (log 10)"))
 csum <- with(df5raw, cumsum(nReads * nOccurrences))</pre>
 csum <- csum / csum[length(csum)]</pre>
 print(xyplot(csum ~log10(nOccurrences), df5raw,
               xlab="Copies per read (log 10)",
               ylab="Cumulative proportion of reads",
               type="1"))
```

Results appear in Figure 2. The cumulative form of this figure appears in the QA report.

The power-law relationship between copies per read and number of unique reads in the control lane consists of three components. At the left of the graph are  $>10^5$  reads that are each represented by only one copy. These likely correspond to sequencing, base calling, or other errors associated with the technology. At the right of the figure are a small number of reads represented many times. These 'frequent' sequences are summarized in the frequentSequences element of qa; frequent sequences are also reported by the tables function of ShortRead.

# Exercise 20

Discover the frequent sequences amongst the raw and aligned reads of lane 5.

Frequent sequences include poly-A reads, reads where only a few bases were called, and reads with close similarity to the Solexa primer or adapter sequence used in sample preparation.

#### Exercise 21

Many of the primer sequences are filtered out by Solexa criteria, but it is worth discovering how many reads are 'similar' to this sequence. Use the srdistance function to identify such reads, e.g., amongst those reads that contain no N nucleotides.

```
> seq <- "CGGTTCAGCAGGAATGCCGAGATCGGAAGACGGT"
> dist <- srdistance(clean(aln), seq)[[1]]
> head(table(dist))

dist
    0    1    2    3    4    5
218 107    59    31    25    17
```

srdistance returns the edit distance between each read and the reference sequence, where the edit distance is defined so that each base mismatch represents an additional increment of 1. There are 384 reads that differ at 2 or fewer locations, from amongst the 370049 reads in the cleaned sample used in this exercise.

Reads represented many times may be problematic for downstream analysis. For instance, sample preparation protocols may involve a PCR step that results in differential amplification, whereas the analysis assumes reads are represented in proportion to their occurrence. The ualignFilter function defined above, and the srduplicated function in *ShortRead* represent two approaches to dealing with this problem by ensuring that reads are represented exactly once (by some definition of 'once'!).

Sequences in the middle portion of the graph in Figure 2 will often, depending on the nature of the investigation, represent the sequences of main biological interest. These are sequences represented an intermediate number of times, as might be required for reasonable coverage in a SNP discovery or ChIP-seq experiment. The right-hand graph in Figure 2 shows a relatively abrupt transition between reads represented rarely and those represented many times.

The sample QA report shows that the non-control lanes show much broader transitions from rarely to frequently represented reads. This could represent technical shortcomings of this run (e.g., inadequate enrichment of sample DNA)

or features of intrinsic biological interest (e.g., wide variability in ChIP abundance between binding sites). Regardless of ultimate source, the broad distribution of read occurrences implies some effort may be required to distinguish 'noise' (reads corresponding to those in the left and right portions of the control lane graph) from signal.

Finally, while the control lane shows a relatively abrupt transition between reads that occur rarely and those that are common (Figure 2), the distribution is in fact 3- or 4-fold broader than expected under a naive model of random read starts along the  $\varphi$ X-174 genome. This is reinforced by alignments to the reference genome, where clear patterns (e.g., unequal representation on plus and minus strands; non-uniform coverage) are apparent.

#### Exercise 22

As an advanced exercise, simulate reads selected uniformly along both strands of the 5200bp long  $\varphi X$ -174 genome. Compare the times each read is represented in your sample with those from the actual control lane. Hint: use sample with replace=TRUE to generate the reads, and table(table(reads)) to summarize their occurrence; this should take less than 5 lines of R code.

# 4.4 Cycle-specific qualities and base calls

As a final foray into the details of quality assessment, consider base calls and quality, and how these change across cycles (see the second table and Section 4 of the QA report).

The table in the QA report suggests that the control lane (lane 5) is enriched for A and T; this is confirmed by the figure in section 4. This likely reflects underlying differences in the genomic regions represented in each lane.

#### Exercise 23

Use alphabetFrequency to summarize nucleotide use in the short reads in aln, from the first part of this lab.

```
> alphabetFrequency(sread(aln), collapse=TRUE, baseOnly=TRUE,
+ as.prob=TRUE)
```

```
A C G T other 0.25588029 0.22092086 0.20225720 0.22800937 0.09293229
```

The frequency of 'other' (i.e., uncalled) nucleotides (> 9%) is very high; typical runs are < 1%.

An unexpected aspect of the figure in Section 4 of the QA report is apparent trends in nucleotide frequency with cycle. For instance, all lanes show a marked decrease in A and increase in C across cycles. This is unexpected in this experiment, where the *a priori* expectation is that sequences start at essentially arbitrary locations in the sequenced DNA: an A is expected as frequently at the beginning of the sequence as at the end.

### Exercise 24

Use alphabetByCycle to extract the number of each nucleotide A, C, T, G, sequenced at each cycle. Plot the transpose of the matrix using matplot, to see how nucleotide counts change across cycles.

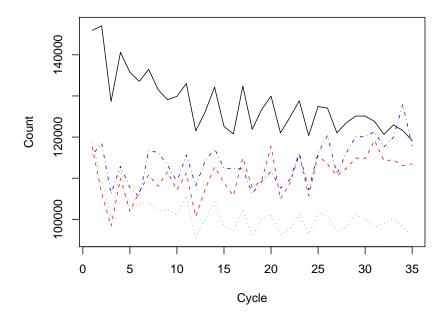


Figure 3: Nucleotide frequency per cycle, lane 1

```
> abc <- alphabetByCycle(sread(aln))</pre>
> dim(abc)
[1] 17 35
> abc[1:4,1:5]
        cycle
                   [,2]
                          [,3]
                                 [,4]
alphabet
           [,1]
                                         [,5]
       A 145955 147044 128640 140609 135751
       C 117538 106357 98356 109948 101961
       G 118619 109879 100303 112267 102233
       T 115666 118470 106223 112917 107716
> matplot(t(abc[1:4,]), type="l", xlab="Cycle", ylab="Count")
```

There are a number of possible contributors to cycle-dependent nucleotide frequencies, including inadequate reagent volume, and nucleotide-specific differential accumulation of fluorescent dyes. Figure 3 and the figure in Section 4 of the report contain additional features that are moderately unexpected, and unexplained. For instance, the frequency of a nucleotide such as A changes very systematically across cycles, first increasing and then decreasing; this seems more regular than expected. Patterns of nucleotide change also seem to echo one another at similar cycles but in different lanes, even when the lanes have different biological material. This occurs for instance in lanes 1-4 of the QA

report, where the last 5 cycles of the C nucleotide seem to change in (comparative!) unison.

# 5 Coverage

Many work flows invove some aspect of coverage, the depth to which individual nucleotides in a reference sequence are covered by nucleotides in reads.

#### Exercise 25

Select just those reads that aligned to the reference genome (e.g., using !is.na(position(aln))). Use the coverage function to summarize how reads are aligned to chromosomes

```
> cvg <- coverage(aln[!is.na(position(aln))])</pre>
```

The coverage function returns a list-like structure, where each element of the list represents a chromosome. The data are represented as a 'run length encoded' (Rle) instance. Rle instances are readily interogatted for a variety of useful insights.

### Exercise 26

Select an element of cvg for further investigation. Use (slice) to identify all regions where there are at least 2 reads. Use the width accessor and density-plot function to display the distribution of island widths. Use viewMaxs and desnityplot to display the distribution of maximum island depths.

```
> slc <- slice(cvg[["chr1.fa"]], 2)
> densityplot(log10(width(slc)))
> densityplot(viewMaxs(slc))
```

There are many additional exploratory opportunities associated with coverage vectors; we will encounter some of these in ChIP-seq processing. As a final exercise. . .

## Exercise 27

The rtracklayer package can export Rle and other objects to a format such as gff suitable for viewing in genome browsers. Use rtracklayer and its export function to export a description of the peaks.

```
> library(rtracklayer)
> fl <- tempfile()
> export(cvg[["chr1.fa"]], fl, "gff3")
```

# 6 Summary

This portion of the course has provided you with an overview of the input, manipulation, and quality assessment functionality available in the *ShortRead* and other packages. A summary of the insights learned might be reflected in simple, and somewhat naive, work flows.

The first work flow simply performs quality assessment on ELAND aligned data.

```
> rpt <- report(qa(extdataDir, "_export.txt"),
+ dest="reports/my_report.pdf")</pre>
```

Performing QA on ELAND data does not commit us to using ELAND alignments in subsequent steps.

The second work flow reads aligned data in to R. The work flow might start with aligned reads created by one of many aligners; we start with ELAND \_export.txt files. There are many possible issues highlighted in the forgoing discussion. We choose to establish a series of filters to eliminate some reads at the very start of our work flow. We eliminate reads with ambiguous base calls and failing Solexa's internal filtering criteria. Reads aligning to multiple locations in the genome are not straight-forward to deal with, and are not essential for ChIP-seq style experiments (this is not necessarily the case for expression or RNA-seq experiments), so we remove these. Most close matches to the Solexa primer sequence are flagged as not passing Solexa base call filters, but we eliminate reads near to this as well. Finally, we restrict our attention to those reads that align to assembled nuclear chromosomes, putting aside for the moment those reads aligning to organelle genomes (the X and Y chromosomes also require special consideration, and we might often eliminate these from a first-pass work flow, too). Our second work flow is thus:

```
> filt1 <- nFilter()
> filt2 <- alignDataFilter(expression(filtering=="Y"))
> filt3 <- alignQualityFilter(threshold=1)
> filt4 <- srdistanceFilter("CGGTTCAGCAGGAATGCCGAGATCGGAAGAGCGGT", 4)
> filt5 <- chromosomeFilter("chr[0-9]+.fa")
> filt <- compose(filt1, filt2, filt3, filt4, filt5)
> aln <- readAligned(extdataDir, pattern, filter=filt)</pre>
```

This work flow applies equally to MAQ aligned data, with the exception that Solexa filtering criteria are not available and the chromosome naming convention in the MAQ-aligned reads in our sample are different:

```
> maqDir <- "path/to/maq"
> filt5 <- chromosomeFilter("chr[0-9XY]+$")
> filt <- compose(filt1, filt3, filt4, filt5)
> maq <- readAligned(maqDir, "s_8.map", "MAQMap", filter=filt)</pre>
```

Each of the filters represents a decision. The decision may be inappropriate for particular analyses, and may be revisited as understanding of the data matures.

The *GappedAlignments* class and the *GenomicRanges* and *Rsamtools* packages provide comparable facilities for input and manipulation of reads aligned with indels.

# 7 Session information

- R version 2.11.1 Patched (2010-05-31 r52167), i386-apple-darwin9.8.0
- Locale: C/C/C/C/en\_US.UTF-8
- Base packages: base, datasets, gr<br/>Devices, graphics, methods, stats, tools, utils

- Other packages: AnnotationDbi 1.10.1, BSgenome 1.16.5, BSgenome.Hsapiens.UCSC.hg19 1.3.16, BSgenome.Scerevisiae.UCSC.sacCer2 1.3.16, Biobase 2.8.0, Biostrings 2.16.6, CSAMA10 0.0.3, EatonEtAlChIPseq 0.0.1, GenomicFeatures 1.0.3, GenomicRanges 1.0.5, IRanges 1.6.8, RCurl 1.4-2, Rsamtools 1.0.5, SNPlocs.Hsapiens.dbSNP.20090506 0.99.1, ShortRead 1.6.2, biomaRt 2.4.0, bitops 1.0-4.1, chipseq 0.4.0, hgu95av2probe 2.6.0, lattice 0.18-8, rtracklayer 1.8.1
- Loaded via a name space (and not attached): DBI 0.2-5, RSQLite 0.9-1, XML 3.1-0, grid 2.11.1, hwriter 1.2