

Metabolomics data analysis with Bioconductor

Johannes Rainer (Eurac Research, Italy)¹

June 12, 2017 @CSAMA2017

¹email: johannes.rainer@eurac.edu, github/twitter: [jotsetung](#) 

Talk content

- Very short introduction to metabolomics data analysis.
- Focus on pre-processing of LCMS data.
- Focus on the `xcms` package (*new* user interface), but other exist too (e.g. `yamss`).


Metabolomics?

- Is the large-scale study of small molecules (metabolites) in a system (cell, tissue or organism).
- Metabolites are intermediates and products of cellular processes (metabolism).
- Metabolome?:
 - **Genome**: what can happen.
 - **Transcriptome**: what appears to be happening.
 - **Proteome**: what makes it happen.
 - **Metabolome**: what actually happened. Influenced by genetic and environmental factors.

How are we measuring that?


- Nuclear magnetic Resonance (NMR) - not covered here.
- Mass spec (MS)-based metabolomics

Mass Spectrometry (MS)

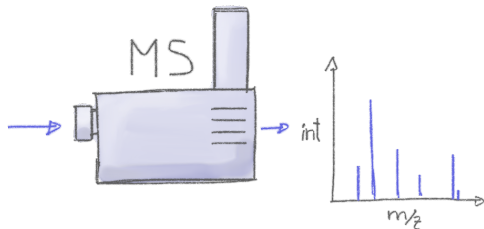

- Problem: unable to distinguish between metabolites with the same mass-to-charge ratio (m/z).

Liquid Chromatography Mass Spectrometry (LCMS)

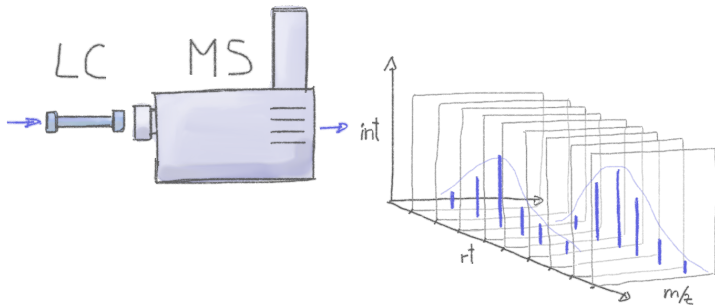

- Combines physical separation via LC with MS for mass analysis.
- Additional time dimension to separate different ions with same m/z .
- LCMS metabolomics: identify peaks in the m/z - rt plane.

LCMS-based metabolomics data pre-processing

- **Input:** mzML or netCDF files with multiple MS spectra per sample.
- **Output:** matrix of abundances, rows being *features*, columns samples.
- **feature:** ion with a unique mass-to-charge ratio (m/z) and retention time.
- **Example:** load files from the faahK0 data packages, process using xcms.

```
library(xcms)
library(faahK0)
library(RColorBrewer)

cdf_files <- dir(system.file("cdf", package = "faahK0"), recursive = TRUE,
                full.names = TRUE)[c(1, 2, 7, 8)]

## Read the data
faahK0 <- readMSData2(cdf_files)
```

- OnDiskMSnExp: small memory size, loads data on-demand.

LCMS-based metabolomics data pre-processing


- Chromatographic peak detection.
- Sample alignment.
- Correspondence.

LCMS pre-processing: Peak detection

- **Goal:** Identify chromatographic peaks within slices along mz dimension.
- What type of peaks have to be detected?


```
mzr <- c(241.1, 241.2)
chr <- extractChromatograms(faahK0, mz = mzr, rt = c(3550, 3800))

cols <- brewer.pal(3, "Set1")[c(1, 1, 2, 2)]
plotChromatogram(chr, col = paste0(cols, 80))
```

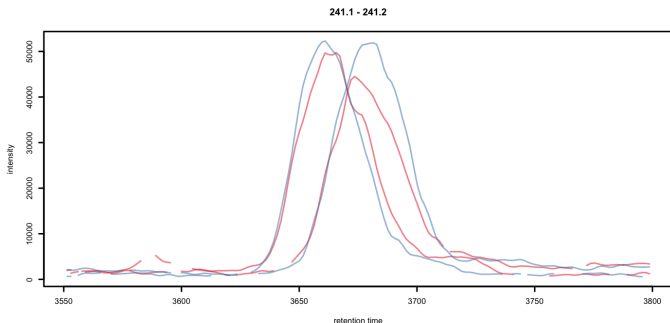

LCMS pre-processing: Peak detection


- **centWave** (Tautenhahn et al. *BMC Bioinformatics*, 2008):
- Step 1: Detection of regions of interest

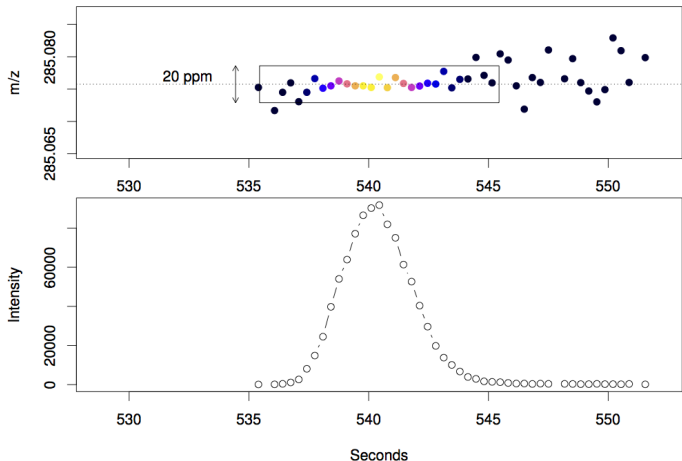

- **mz-rt regions with low mz-variance.**

LCMS pre-processing: Peak detection

- Step 2: Peak detection using continuous wavelet transform (CWT)

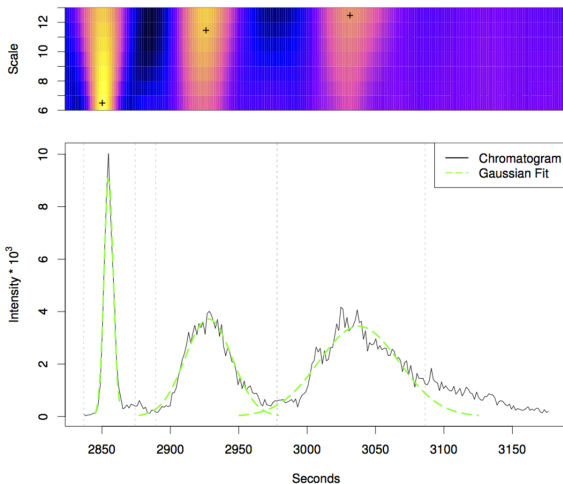

- Allows to identify peaks with different widths.

LCMS pre-processing: Peak detection

- Example: centWave-based peak detection:

```
faahK0 <- findChromPeaks(faahK0, param = CentWaveParam())
```

- Result: XCMSnExp, container for LC/GC-MS results, extends OnDiskMSnExp.


```
head(chromPeaks(faahK0))
```

	mz	mzmin	mzmax	rt	rtmin	rtmax	into	intb	maxo	
[1,]	425.9	425.9	425.9	425.9	2520.158	2510.768	2527.982	9999.769	9984.120	741
[2,]	464.3	464.3	464.3	464.3	2518.593	2504.508	2532.677	32103.270	32082.926	1993
[3,]	499.1	499.1	499.1	499.1	2524.852	2520.158	2527.982	4979.194	4904.709	883
[4,]	572.7	572.7	572.7	572.7	2524.852	2520.158	2527.982	2727.446	2721.187	559
[5,]	579.8	579.8	579.8	579.8	2524.852	2520.158	2527.982	2450.477	2444.218	468
[6,]	453.2	453.2	453.2	453.2	2506.073	2501.378	2527.982	1007408.973	1007380.804	38152

	sn	sample	is_filled
[1,]	740	1	0
[2,]	1992	1	0
[3,]	13	1	0
[4,]	558	1	0
[5,]	467	1	0
[6,]	38151	1	0

LCMS pre-processing: Alignment

- **Goal:** Adjust retention time differences/shifts between samples.
- Total ion chromatogram (TIC) representing the sum of intensities across a spectrum.

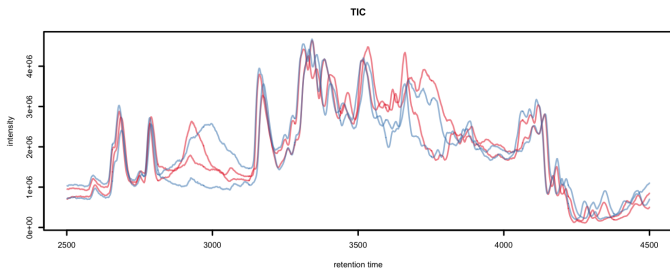


- Overview of algorithms: (Smith et al. *Brief Bioinformatics* 2013).
- *xcms*: *peak groups* (Smith et. al *Anal Chem* 2006), *obiwarp* (Prince et al. *Anal Chem*, 2006),

LCMS pre-processing: Alignment

- Example: use obiwrap to align samples.

```
faahK0 <- adjustRtime(faahK0, param = ObiwrapParam())
```

- TIC after adjustment:

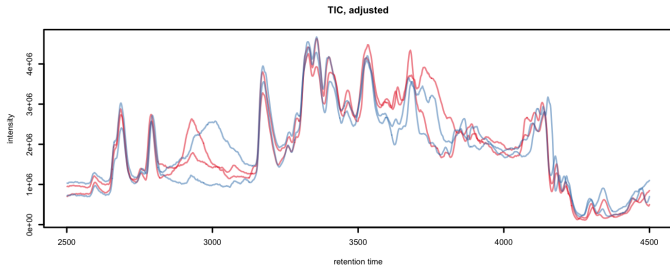

- Assumptions:
 - Samples relatively similar (either similar chromatograms or a set of common metabolites present in all).
 - Analyte elution order same in all samples.

LCMS pre-processing: Alignment



- Example: effect of alignment on example peak.

```
chrs_adj <- extractChromatograms(faahKO, mz = mzr, rt = c(3550, 3800))
```

```
par(mfrow = c(2, 1))
```


```
plotChromatogram(chrs, col = paste0(cols, 80), main = "Before alignment")
```

```
plotChromatogram(chrs_adj, col = paste0(cols, 80), main = "After alignment")
```

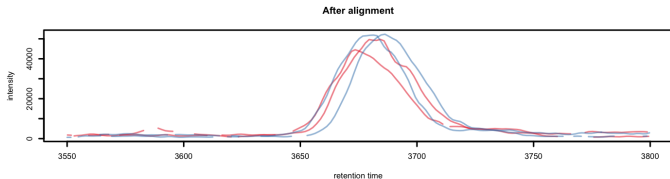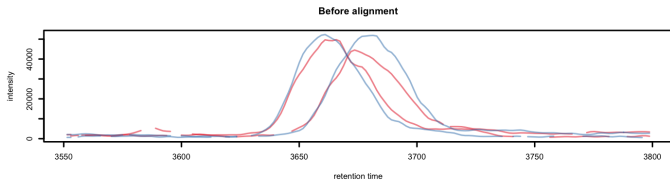

LCMS pre-processing: Correspondence

- **Goal:** Group detected chromatographic peaks across samples.
- Peaks that are close in rt (and m/z) are grouped to a *feature*.
- xcms: *peak density* method:

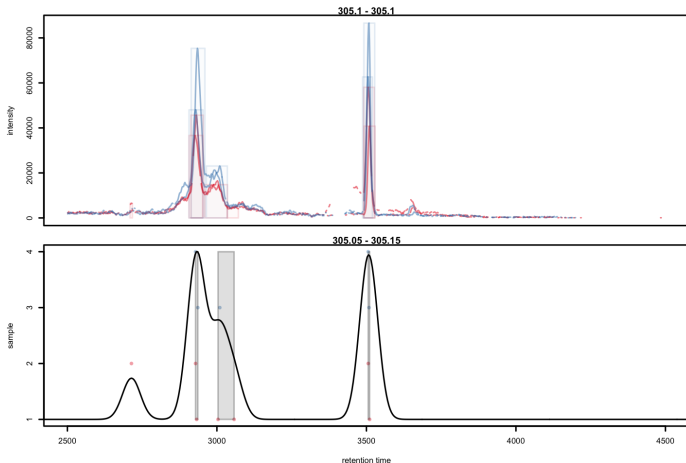

LCMS pre-processing: Correspondence

- Example: peak grouping.

```
faahK0 <- groupChromPeaks(faahK0, param = PeakDensityParam())
```

- `featureValues`: extract *values* for each feature from each sample.


```
## Access feature intensities  
head(featureValues(faahK0, value = "into"))
```

	ko15.CDF	ko16.CDF	wt15.CDF	wt16.CDF
FT0001	6029.945	NA	4586.527	NA
FT0002	1144.015	NA	1018.815	NA
FT0003	NA	774.576	1275.475	NA
FT0004	NA	NA	1284.728	1220.7
FT0005	2759.095	3872.963	NA	NA
FT0006	7682.585	3806.080	NA	NA

- *Fill-in* values for missing peaks: `fillChromPeaks`.

What next? Data normalization

- Adjust within and between batch differences.
- MetNorm RUV for metabolomics (Livera et al. *Anal Chem* 2015).
- Injection order dependent signal drift (Wehrens et al. *Metabolomics* 2016).

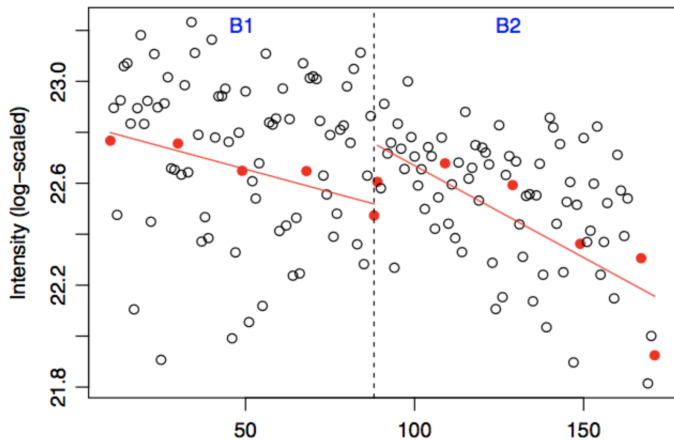

What next? Identification

- Annotate features to metabolites.
- Each metabolite can be represented by multiple features (ion adducts, isotopes).
- Starting point: CAMERA package.
- On-line spectra databases (e.g. MassBank).

Finally. . .

thank you for your attention!

- Hands on in the afternoon labs:
 - Proteomics lab.
 - Metabolomics lab (pre-processing of LCMS data).